Hot or Cold Pack

REGULAR MODEL EN

PRODUCT FEATURES

Three-dimensional structure, more temperature control materials.

Longer cooling/heating time

No risk of leakage (210D TPU)

No tilting of the inner material to one side during use

HOW TO USE

- **COLD PACK:** Place in the freezer 0°F (-18°C) for about 3 hours, Then remove and apply to the desired cold compress area.
- **HOT PACK:** Heat in the microwave 122°F (50°C) for 30 seconds (The number of heating times can be adjusted according to the desired temperature). Once the temperature is suitable, apply to the affected area for relief.r relief.

TIME OF USE

COMPARISON WITH PEERS			
ENVIRONMENTAL TEMPERATURE (77°F/25°C)			
Temperature	≤50°F/10°C	≤59°F/15°C	≤68°F/20°C
Our	249mins	274mins	318mins
Other	50mins	82mins	151mins
Temperature	32-50°F/(0-10°C)	32-59°F/(0-15°C)	32-68°F(0-20°C)
Our	216mins	242mins	286mins
Other	28mins	60mins	129mins

Test temperature: 77°F(25°C) constant temperature;

Test tools: constant temperature and humidity chamber, digital thermometer, hot or cold pack, 2mm pearl cotton;

Thermometer probe position: hot or cold pack is placed on pearl cotton, and the thermometer probe is placed between hot or cold pack and pearl cotton.

Data display: Data position 77°F(25°C) ambient temperature, hot or cold pack duration in different temperature ranges.

Note: The test is a simulation test, not a real human test.

WORKING PRINCIPLE

- Principle of cold compress: By lowering the surface temperature of the skin and slowing down the blood flow, it plays the role of anti-inflammatory, pain relief and reducing swelling.
- Principle of Hot Compress: By local warming, it promotes blood circulation, accelerates tissue repair, and relieves pain and muscle tension.

USAGE SCENARIOS

- Daily body therapy.
- Cold therapy for muscle strains and joint pain after exercise.
- Hot compress therapy for chronic pain (e.g., herniated discs, arthritis).
- Heat therapy to relieve abdominal pain during menstruation.
- Emergency cold compress for everyday bruises and injuries.
- Relief of shoulder and neck pain after long hours of desk work.
- Muscle relaxation and recovery after exercise.

CAUTIONS

- Before use, please check the product for any damage and ensure there are no leaks.
- After cold therapy, avoid engaging in intense physical activity immediately.
- When using for hot therapy, be mindful of the temperature to avoid burns.
- Children should use the product under adult supervision.

WARNINGS

- Do not place the product directly on an open flame or high-temperature heaters.
- Do not puncture the product or ingest its contents.
- Do not continue to use the product if it is damaged. Avoid contact between the contents and skin or eyes.

CONTRAINDICATIONS

- Not for use on areas with open wounds or skin infections.
- Strictly prohibited for uninstructed use in people with loss of temperature sensitivity (e.g. diabetics).
- Pregnant women should consult a physician before use.

FAQS

Q: Is the product washable?

A: You can use a damp cloth to gently wipe the surface, avoid direct washing.

Q: What should I do if the product is not effective in cold application?

A: Make sure the product is well chilled in the freezer 0°F (-18°C) for at least 3 hours.

Q: What if the temperature is too high after microwave heating?

A: It is recommended to test the temperature after each heating, and if necessary, wait until it cools down to a suitable temperature before use.

PRODUCT INFORMATION

Brand: TEMPK™

Name: Hot or Cold Pack - Regular Model

Size: 255*130*20mm (±5)

Material of outer bag: 210D TPU

Contents: CMC, Water, Glycerol, Preservative

Safety: The products are made of environmentally friendly materials, non-toxic and harmless

Producer: Shanghai Huizhou Industrial Co., Ltd.


Origin: China

PACKAGE CONTAINS

x 2

Hot/Cold Pack

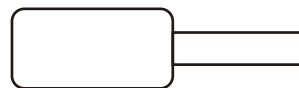x 1

Velvet Cover

CONTACT US

If you need assistance, please contact us.

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Fabricante: Shanghai Huizhou Industrial Co., Ltd.